

SAFE SHARPS PRODUCTS ORDERING INFORMATION

Please note: This document contains examples of safe sharps devices and is for informational purposes only. EHS does not endorse any of these products, but is recommending consideration of their use for safety reasons. Other companies may have similar products.

Safety Needles**General Tips:**

- Immediately dispose needle and/or syringe after use in an approved sharps container.
- Always keep a sharps container in immediate vicinity of use for direct disposal.
- Avoid recapping, bending or breaking needles.

VanishPoint Syringe
with retractable needle,
Retractable Tech,
various sizes, for 23G x
1 inch

SciQuest (McKesson) Part
Number: **348827BX**

Manufacturer (Retractable
Tech) Part Number: **10311**




[VanishPoint Retractable Safety Needle](#)

- Retractable needle
- Different sizes and volumes available

BD Integra syringe with
retractable needle,
various sizes, for 25G x
1 inch

SciQuest Part Number (Fisher
Scientific): **1482351**

Manufacturer (BD Vacutainer
Labware Medical) Part Number:
305270

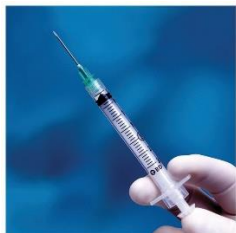

[BD Integra General-Use Syringes](#)

- Retractable needle
- Different sizes and volumes available

BD Luer-Lok Syringe with Detachable Needle 25G x 0.625 inch SciQuest Part Number: BD305781 Manufacturer (Avantor) Part Number: 305781		BD Luer-Lok Syringe with Detachable Needle	• Pivoting shield that covers needle point after use • Luer-lock connection • Different sizes and volumes available
--	---	--	---

Safe Scalpels

General Tips:

- Always keep the blade facing away from you (even if you do not hold it).
- Always keep your hands out of the path of the scalpel.
- Use cut resistant gloves when working with scalpels.
- Never remove scalpel blades by hand.

Excelta Thumb Scalpel Blade Size 11 SciQuest Part Number: 94024-138 Manufacturer (VWR) Part Number: 175-11		Excelta Corp Scalpel with Blade	• Mini scalpel for fine dissection procedures • Single-use • Direct disposal
Disposable Mini-Scalpels Blade Size 10 SciQuest Part Number: N/A Manufacturer (AGRIsales inc.) Part Number: V18362245		Disposable Mini-Scalpels - #10	• Mini scalpel for fine dissection procedures • Single-use • Direct disposal

DeRoyal Retractable Safety Scalpel Swann-Morton Retractable Safety Scalpel Blade Size 11 SciQuest Part Number: QTXD4511CSZ Manufacturer (Medline) Part Number: D4511		DeRoyal Retractable Safety Scalpels	• Retractable blade • Disposable scalpel • Different blade sizes available (color-coded)
Aspen Surgical Bard-Parker Protected Disposable Scalpel Blade Size 10 SciQuest (Fisher Scientific) Part Number: 0268878 Manufacturer (Aspen Surgical) Part Number: 372610		Aspen Surgical Bard-Parker Protected Disposable Scalpel	• Protected blade/Retractable cover • Disposable scalpel • Different blade sizes and handle lengths available

Safe Scalpels			
Aspen Surgical Bard-Parker Protected Blade Size 10 SciQuest Part Number: B-D373910H Manufacturer (Aspen Surgical) Part Number: 373910		Bard-Parker Safety Blades	• Protected blade/Retractable cover • Disposable blade • Quick release • Different blade sizes available

Aspen Surgical Bard-Parker Scalpel Handle Size 3 SciQuest Part Number: B-D374030 Manufacturer (Aspen Surgical) Part Number: 374030		Bard-Parker Safety Handles	• Reusable handle • Metal handle • Quick release of blade • Different handle lengths available
---	---	--	---

Safe Scalpel Removal Tools

General Tips:

- Always remove blades with an engineering control as seen here (never with your hands).
- Always keep a blade removal tool in vicinity of use for direct disposal.

Single-Use Scalpel Blade Remover SciQuest Part Number: N/A Manufacturer (QLICKSmart) Part Number: QSBS01		QLICKSmart Scalpel Blade Remover	• Single-use • Direct disposal
Aspen Surgical Bard-Parker Blade Remover SciQuest Part Number: 22027983 Manufacturer (Aspen Surgical) Part Number: ST-115		Bard-Parker® Blade Remover Aspen Surgical (aspensurgical.com)	• Single-use (holds up to 5 blades) • Disposable

Safe Scalpel Removal Tools

Scalpel Blade Remover SciQuest (VWR/Avantor) Part Number: 21900-000 Manufacturer (GF Health Products) Part Number: 2990		Scalpel Blade Remover	• Multiple use (holds up to 100 blades) • Disposable
Qlicksmart BladeFLASK Scalpel Blade Remover SciQuest (Fisher Scientific) Part Number: 22444275 Manufacturer (Mopec Inc.) Part Number: AJ150		Mopec Blade FLASK Blade Remover, Qlicksmart	• Multiple use (holds up to 100 blades) • Disposable
Integra Miltex Blade-Safe Remover SciQuest (Fisher Scientific) Part Number: 12460439 Manufacturer (Integra Lifesciences Sales): 4-100		Integra™ Miltex™ Blade-Safe Surgical Blade Remover	• Reusable • Autoclavable • Blade needs to be manually removed from device; always use an engineering control (such as hemostat or forceps) to remove the blade
StatLab Aunex Blade Remover SciQuest (Fisher Scientific) Part Number: STL48CSEA Manufacturer (StatLab) Part Number: 48CS-EA		StatLab Aunex Blade Remover	• Single-handed blade removal • No fingers, forceps, or re-sheathing • Remove and dispose of up to 200 blades

Safe Razor Blade Tools**General Tips:**

- When possible, use single sided (instead of double sided) razor blades to minimize risk of injury
- Always remove blades with an engineering control as seen here (never with your hands)

Stanley Mini-Razor Blade Scraper

SciQuest Part Number: N/A

Manufacturer (Grainger) Part Number: **3Q094, 28-100**

[STANLEY Scraper: Bld. Mat. Carbon Steel, Hdl. Mat. Plastic, Stiff, Straight Bld. Shape](#)

- Reusable
- Razor blade needs to be manually turned; Once the razor blade has been inserted, leave the blade in place and store in a plastic sealed container

Stanley Retractable Razor Blade ScraperSciQuest (Avantor/VWR) Part Number: **300003-958**Manufacturer (ORS NASCO Inc.) Part Number: **680-28-500**

[Stanley Retractable Razor Blade Scraper](#)

- Retractable razor blade
- Reusable

FOSHIO Plastic Razor Blade Scraper

SciQuest Part Number: N/A

Manufacturer (FOSHIO) Part Number: **FM052T2-AU**

[FOSHIO Plastic Razor Blade Scraper Include 2PCS Scraper Tool and 100PCS Blades for Gasket Remover, Labels Decal and Adhesive Remover for Windows and Glass,Blue](#)

- Non-retractable razor blade

Cut Resistant Gloves**General Tips:**

- Purchase cut resistant gloves 1 size smaller than your disposable nitrile gloves.
- Always remember cut resistant gloves are not cut proof.
- Cut resistant gloves will not prevent a needle stick or direct puncture by a sharp object. They are designed to minimize injuries from slicing, scraping, or cutting motions with sharp items. The higher the ANSI cut level rating, the more cut resistant the glove.
- ANSI Standard 105-2016 (test ASTM F2992-15):
 - A6 > 3000 g needed to cut through material with 20 mm of blade travel.
 - A7 > 4000 g needed to cut through material with 20 mm of blade travel.

Ansell Polar Bear Lightweight Dyneema HyFlex Glove, Medium SciQuest (VWR/Avantor) Part Number: 47745-060 Manufacturer (Ansell) Part Number: 222163		HyFlex® 74-045 Lightweight Dyneema® Gloves, White, Ansell	• EHS recommends ANSI cut level 4 gloves as the minimal cut resistant glove for use • Thinner glove liner, offers good dexterity • Reusable
Ansell Polar Bear HyFlex Cut Resistant Glove, Medium SciQuest (Fisher Scientific) Part Number: 19160390C Manufacturer (Ansell) Part Number: 245313		Ansell HyFlex™ 74-048 Lightweight Cut Protection Gloves	• Increased cut-resistance protection (ANSI cut level A6), but dexterity impacted for tasks that require fine mobility skills • Thicker glove liner • Reusable
Hubert Essentials Pro Max Red Dyneema Serrated Cut Resistant Glove, Medium SciQuest Part Number: N/A Manufacturer (Hubert) Part Number: 15571		HUBERT® Essentials Pro Max Red Dyneema Serrated Cut Resistant Glove - Medium	• Increased cut-resistance protection (ANSI cut level A7), but dexterity impacted for tasks that require fine mobility skills • Thicker glove liner • Reusable

Steel Mesh Glove SciQuest Part Number: N/A Manufacturer (ULINE) Part Number: S-18009		Steel Mesh Glove	• Highest cut-resistance protection (ANSI cut level A9), but dexterity impacted for tasks that require fine mobility skills • Thicker glove liner • Reusable
--	---	---	--

Puncture Resistant Gloves

HEXARMOR Knit Cut-Resistant Gloves, Medium SciQuest Part Number: N/A Manufacturer (Grainger) Part Number: 61CU91, 2068X-M		HEXARMOR Knit Gloves	• ANSI cut level A2 on polyurethane palm coating • ANSI-ISEA Abrasion Level 6 • Polyurethane coating provides grip and abrasion resistance • Thinner glove liner, offers good dexterity • Reusable
---	---	---	--

Surgical Tools

Shandon Straight Blunt/Sharp Dissecting Scissors (6.5 in) SciQuest (Fisher Scientific) Part Number: 28251 Manufacturer (Thermo Scientific) Part Number: 28251		Epredia™ Shandon™ Straight Blunt/Sharp Dissecting Scissors, premium, 6.5 in. (16.5cm)	• Safer option for animal necropsy than sharp ended scissors • One tip is blunt and one is sharp
--	---	--	---

Operating Scissors, stainless steel with straight, blunt tips SciQuest (VWR/Avantor) Part Number: 89259-980 Manufacturer (VWR International) Part Number: 89259-980		VWR® Operating Scissors	• Safer option for animal necropsy than sharp-ended scissors. Both tips are blunt
Broad-Tipped Forceps SciQuest (Fisher Scientific) Part Number: 10300 Manufacturer (Fisher Scientific) Part Number: 10300		Fisherbrand™ General-Purpose Broad-Tipped Forceps https://www.fishersci.com/shop/products/fisherbrand-general-purpose-broad-tipped-forceps/10300	• Safer alternative to fine-tipped (sharp) forceps
Eisco Crucible Tongs with Bow, Metal, Length 24cm SciQuest (Fisher Scientific) Part Number: S27100 Manufacturer (Eisco Scientific) Part Number: CH0729B		Eisco™ Crucible Tongs with Bow	• Stainless steel crucible tongs with looped handles • Bowed jaws with flat ends • Can safely grasp hot crucibles, beakers or other apparatus measuring up to 4" in diameter • Tongs measure 8.25" in total length

Surgical Tools			
Baker Original Hookout Hook Remover SciQuest Part Number: N/A Manufacturer (Baker) Part Number: 044434001004		Baker Original Hookout Hook Remover	• Device to hold sharp items for manipulation of needles and syringes when needed (to avoid hand injury) • Reusable

Infinity Reusable OR Towels SciQuest (Medline Industries) Part Number: MDTST5H31CEI Manufacturer (Medline) Part Number: MDTST5H31CEI		Infinity Reusable OR Towels Medline	• Cover sharp items (such as bones) during surgery or autopsy • Reusable/ Washable
Baker 8" Pliers Forceps & Line Clipper/Punch Tool Kit SciQuest Part Number: N/A Manufacturer (Baker) Part Number: B8SSPC		Baker 8" Pliers Forceps & Line Clipper/Punch Tool Kit Tackle Warehouse	• Device to hold sharp items for manipulation of needles and syringes when needed (to avoid hand injury) • Reusable
Plastic Tip Tweezers SciQuest Part Number: 23A00A161 Manufacturer (Ideal-tek SA) Part Number: 259COR.SA.1.ITU		Plastic Repl CO Tip Tweezers	• Safer alternative to fine tipped (sharp) forceps • Plastic tip tweezers • High heat capability (230 – 260 °C) • High chemical resistance to hydrochloric acid and oxidizing agents • Reusable
Medique Plastic Forceps SciQuest (Fisher Scientific) Part Number: 19090651 Manufacturer (Medique Products) Part Number: 77233		Medique Plastic Tweezers	• Economy, disposable plastic tweezer, size 4-1/2"

For assistance with Workday/SciQuest, see [Procurement Toolkit, Requisition to Purchase Order](#)